



1. Navigate the ribbon bar to the patient menu. Select the “binoculars” to access the patient advanced search area.

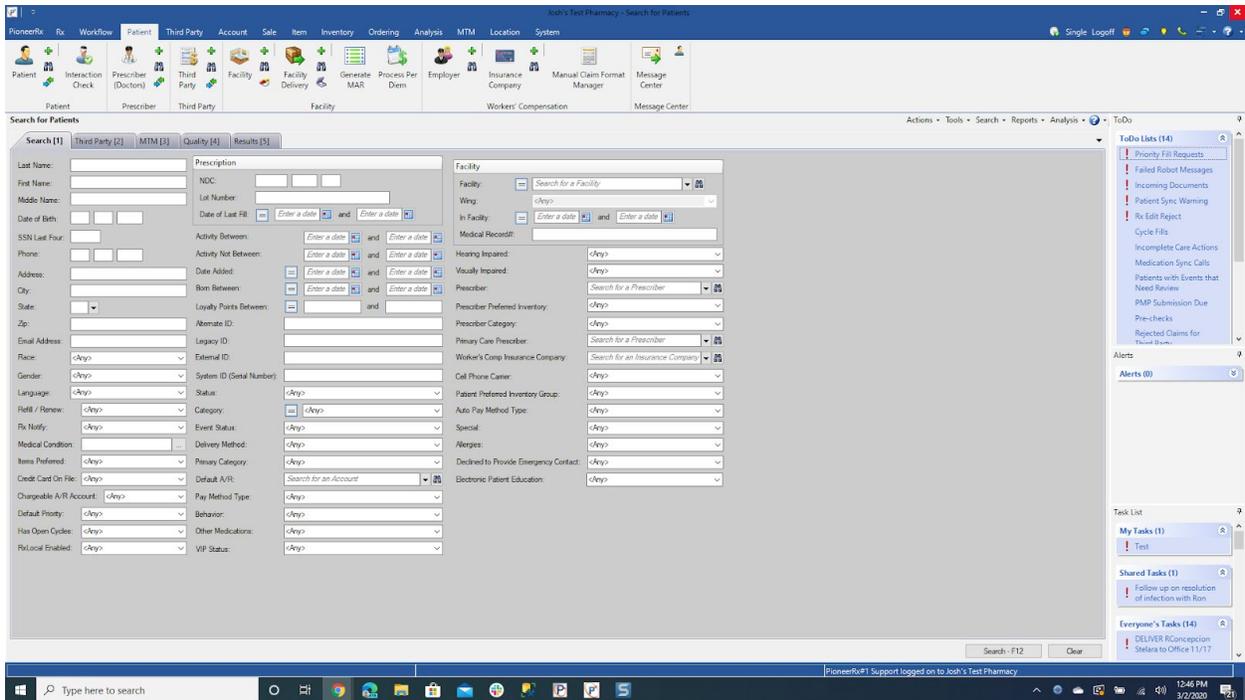

2. Within the patient advanced search area, navigate to the “Reports menu”

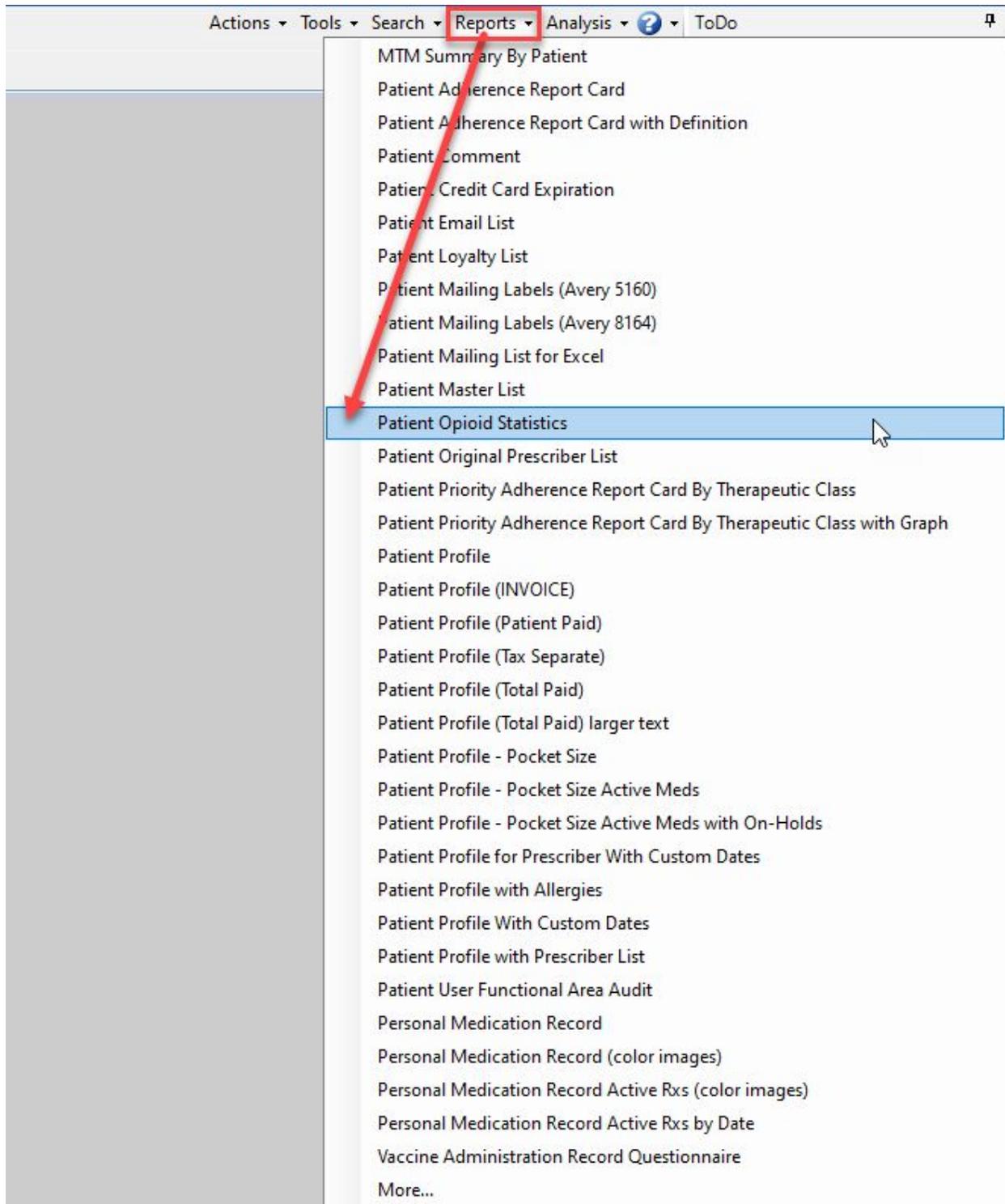

3. Select the "Patient Opioid Statistics" report. Once selected PioneerRx will run for a few minutes while querying the previous 3 month period of time (current date - 90 days). NOTE: if the user enters a date on the patient advanced search area that date will be plugged in instead.

Patient Opioid Statistics Report Printer: EPSON WF-C8690 Series-F15E8A

Print Immediately Print... Fax Excel Text 100 % 1/1 Backward Forward

Patient Opioid Statistics

	Total Unique Patients*	Percentage (Compared to On Opioids)
On Opioids	246	N/A
High Dose Opioids: Concurrently Using Opioids and Benzodiazepine (for >30 Consecutive Days):	11	4.47%
Concurrently Using 2 Short Acting Opioids (for >30 Consecutive Days):	23	9.35%
One or More Opioid Fills in a Consecutive 4-month Period (Chronic Pain):	115	46.75%
	60	24.39%

Total Number of Naloxone Transactions Dispensed Over Last 6 Months: 14

The report has completed

Default report without date filled between and status of active.

Patient Opioid Statistics Report Printer: EPSON WF-C8690 Series-F15E8A

Print Immediately Print... Fax Excel Text 100 % 1/1 Backward Forward

Patient Opioid Statistics

Last Fill Between: 1/1/2020 and Any
Patient Status: Active

	Total Unique Patients	Percentage (Compared to On Opioids)
On Opioids	239	N/A
High Dose Opioids:	8	3.35%
Concurrently Using Opioids and Benzodiazepine (for >30 Consecutive Days):	22	9.21%
Concurrently Using 2 Short Acting Opioids (for >30 Consecutive Days):	115	48.12%
One or More Opioid Fills in a Consecutive 4-month Period (Chronic Pain):	60	25.10%

Total Number of Naloxone Transactions Dispensed Over Last 6 Months: 14

The report has completed

Report with custom date filled

- The summary report is generated. In the area listed “Last fill between 1/1/20 and Any” and “Patient status: Active” are the criteria used on the patient advanced search grid. If no other criteria are entered, the report will run (Current date - 90 days). Current date - 90 days is the default.